


Anderson Pipette

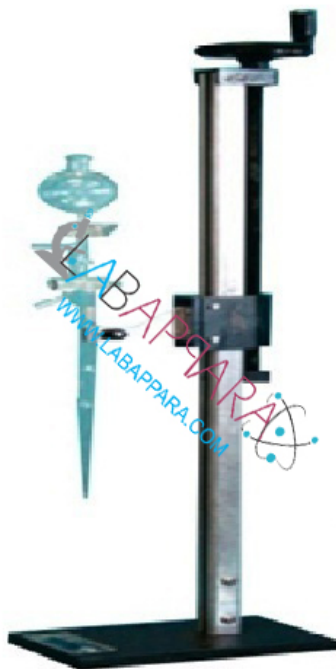

Product Categories: [Laboratory Equipments](#), [Pharmaceutical Equipments](#)

Product Tags: [biology lab equipment](#), [chemistry lab instruments](#)

Product Page: <https://www.labappara.com/product/anderson-pipette/>

Product Description

Anderson Pipette - Laboratory Equipment

Anderson Pipette that is used in laboratories to find out particles size in any suspension.

Particle size Determination apparatus, (Sedimentation pipette, Anderson) consisting of 500ml.

Provided pipette is precisely designed using quality assured raw material and technology in with industry quality standards.

Graduated cylinder, graduation from 0 to 20 cm and fitted with 10 ml.

Pipette with three way stopcock, the tip of pipette at the level of 0 mark at the cylinder.

Product data have been exported from - [Labappara scientific instrument](#)
Export date: Tue Nov 12 16:24:24 2024 / +0000 GMT