


Hamming Code Trainer

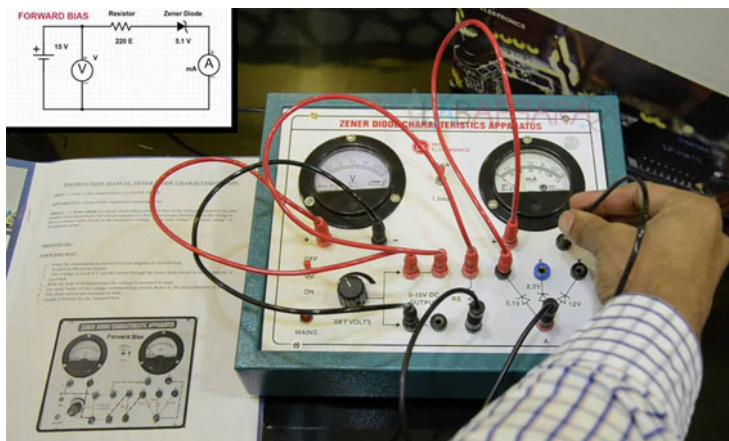

Product Categories: [Advance Communication Lab](#), [Engineering Equipment](#)

Product Page: <https://www.labappara.com/product/hamming-code-trainer/>

Product Description

Hamming Code Trainer

Features:

Rugged design

Light weight

Enhanced performance

Other Details:

Onboard transmitter clock and Sync

(Fast and Slow mode)

On Board Modulating Digital Data signal generator to generate any binary input word with Bit clock & Word Clock.

Word Length: 7 Bits

On Board check bits: Three bits and Four bits data

On Board Noise inserting circuit through four bit error generator

Onboard seven bit shift register

To generate digital data word with four bit data & three parity check bits